






Styleitaliano DAILY by HORICO

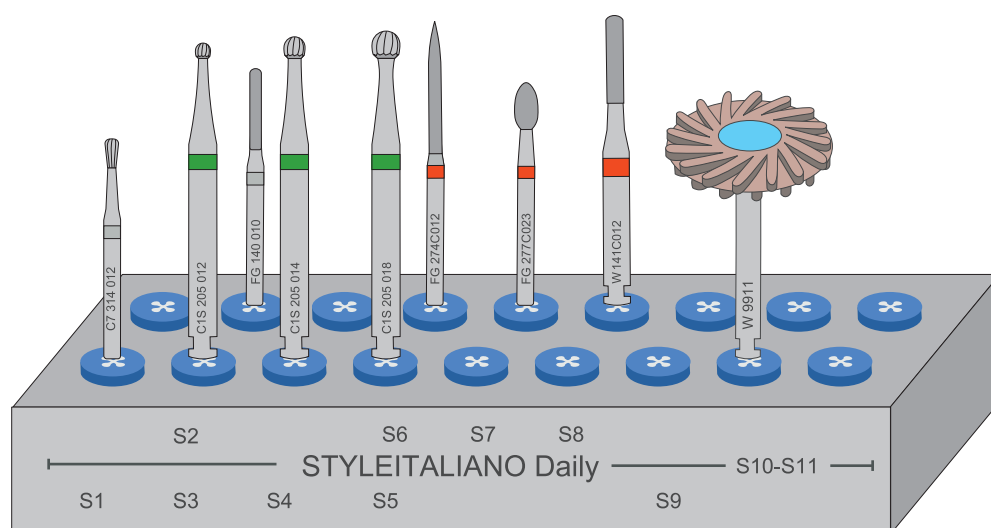





Introduction

The Styleitaliano community, renowned for simplifying and enhancing restorative dentistry, developed with HORICO a specialized bur kit tailored for direct restorations. This comprehensive kit offers a systematic approach with distinct sequences for direct anterior and posterior procedures, ensuring precision and efficiency. The use of protocols makes this kit the best resource for the clinicians and for the dental assistants; easy to organize and reorder.

A standout feature of this kit is its proprietary finishing bur together with the Lamello® polisher is designed to deliver exceptional lifelike surfaces. Complemented by separating diamond strips that enhance the proximal finishing process. This collaboration exemplifies Styleitaliano's commitment to creating tools that are feasible, teachable, and repeatable, streamlining restorative procedures for dental professionals.

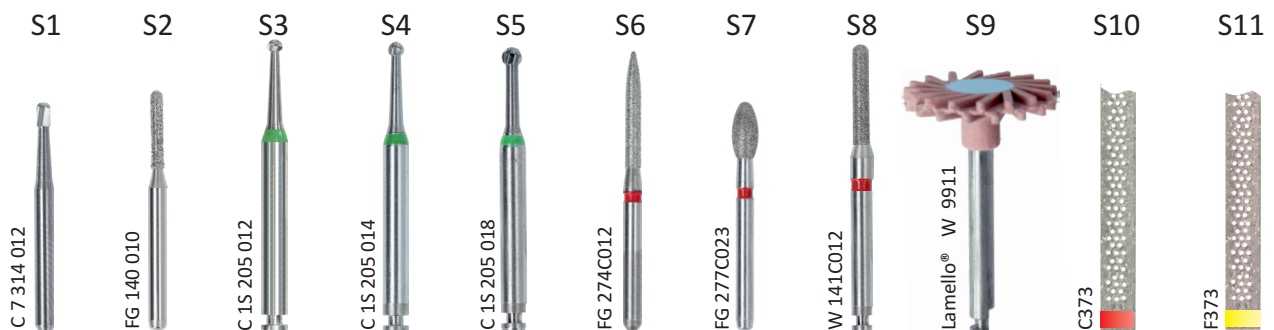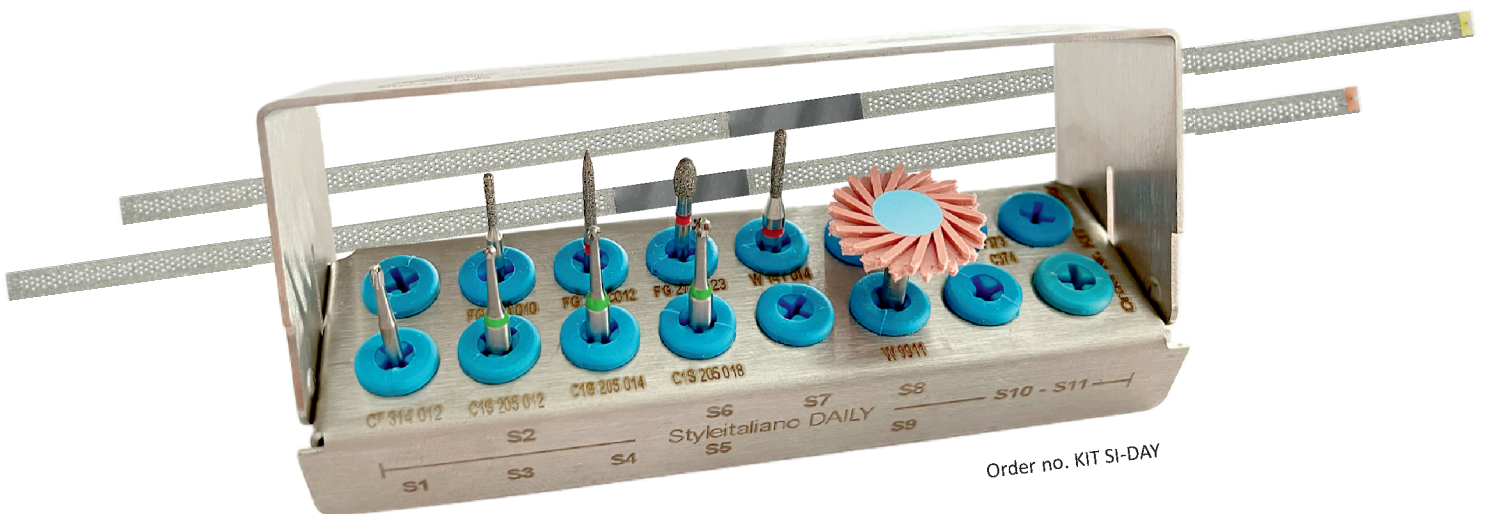


Anterior Sequence

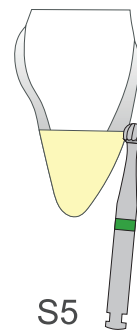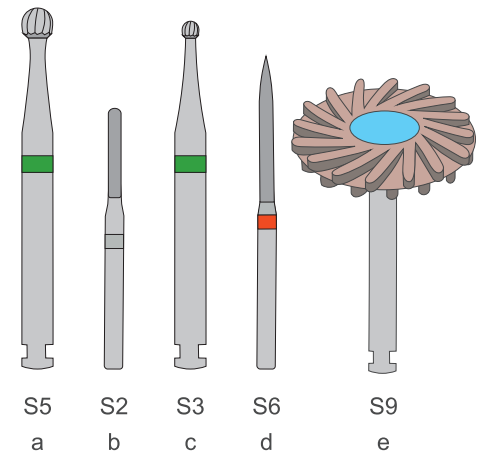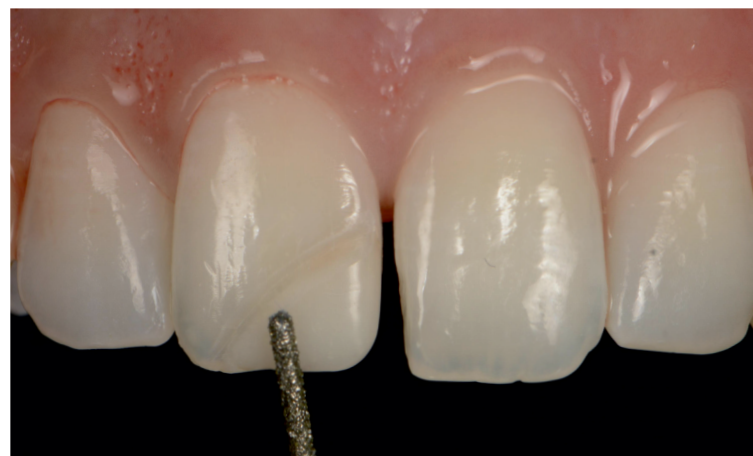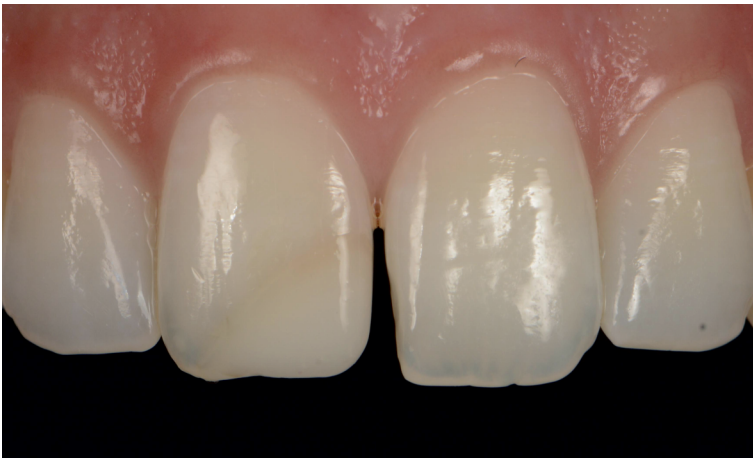

Anterior sequence


- a- Margin finder
b- Restoration removal and prep
c- Deep cleaning
d- Bevel
e- Margin polishing

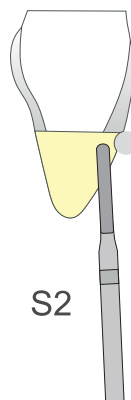

Anterior sequence



- a- Margin finder
b- Restoration removal and prep
c- Deep cleaning
d- Bevel
e- Margin polishing

Anterior Sequence

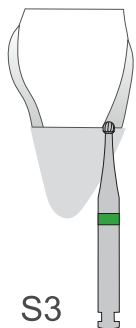


Anterior sequence

- a- Margin finder
- b- Restoration removal and prep
- ➔ c- **Deep cleaning**
- d- Bevel
- e- Margin polishing

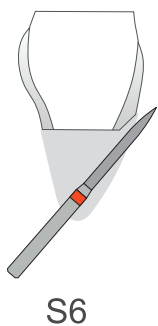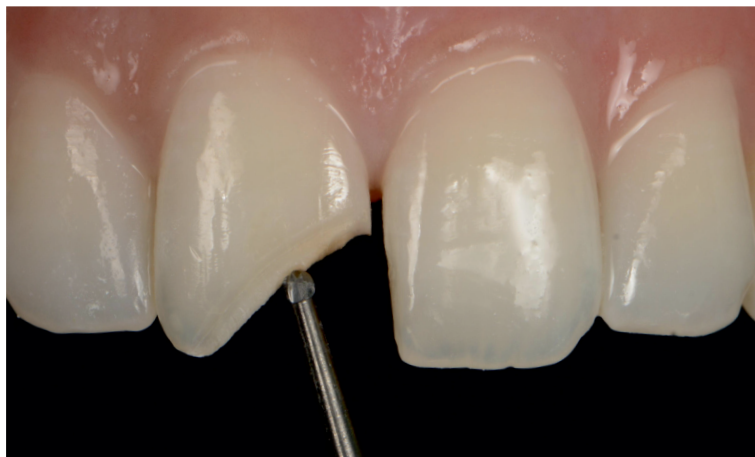

Anterior sequence

- a- Margin finder
- b- Restoration removal and prep
- c- Deep cleaning
- ➔ d- **Bevel**
- e- Margin polishing

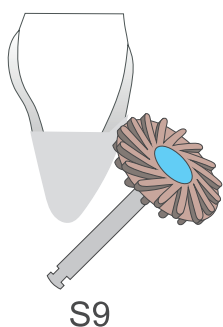

Anterior sequence


- a- Margin finder
- b- Restoration removal and prep
- c- Deep cleaning
- d- Bevel
- ➔ e- **Margin polishing**



Anterior Finishing Sequence




Anterior Finishing Sequence

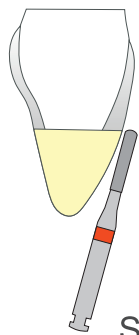

S8

Anterior finishing sequence

- a- Palatal wall correction
- ➔ b- Secondary (and tertiary anatomy)
- c- Pre finishing

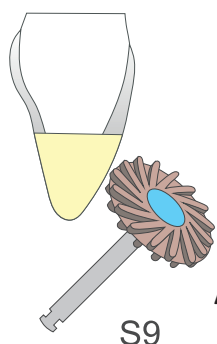

S9

Anterior finishing sequence

- a- Palatal wall correction
- b- Secondary (and tertiary anatomy)
- ➔ c- Pre-finishing




Anterior Finishing Sequence





Anterior Finishing sequence

- a- Palatal wall correction
- b- Secondary anatomy
- c- Pre-finishing

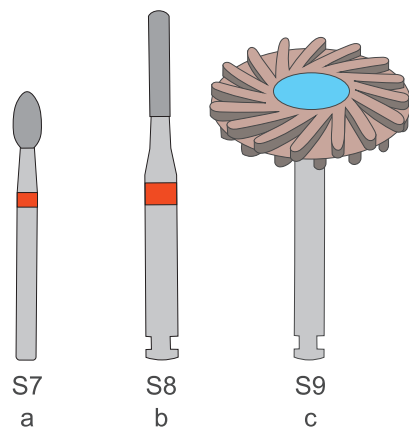


Posterior Sequence

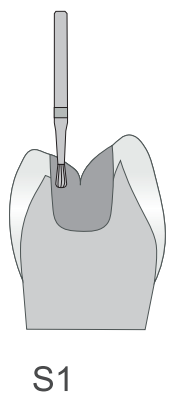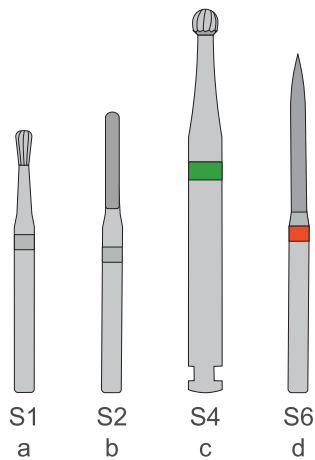

Posterior sequence

- a- Amalgam removal
b- Cavity prep
c- Caries cleaning
d- Box opening and margin smoothing

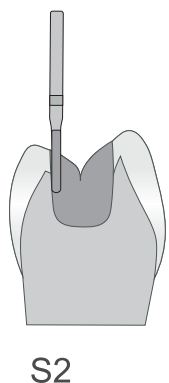


Posterior sequence

- a- Amalgam removal
b- Cavity prep
c- Caries cleaning
d- Box opening and margin smoothing



Posterior Sequence

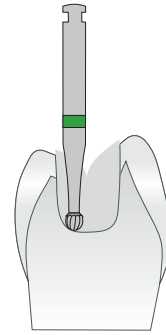

S3 or S4 or S5




Posterior sequence

a- Amalgam removal

b- Cavity prep

c- **Caries cleaning**

d- Box opening and
margin smoothing

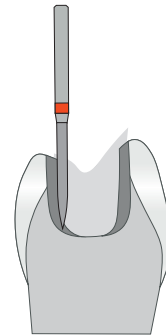

S6



Posterior sequence

a- Amalgam removal






b- Cavity prep

c- Caries cleaning

d- **Box opening and
margin smoothing**



Posterior Finishing Sequence

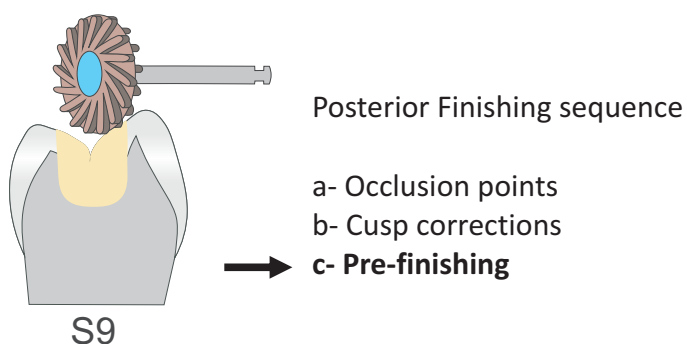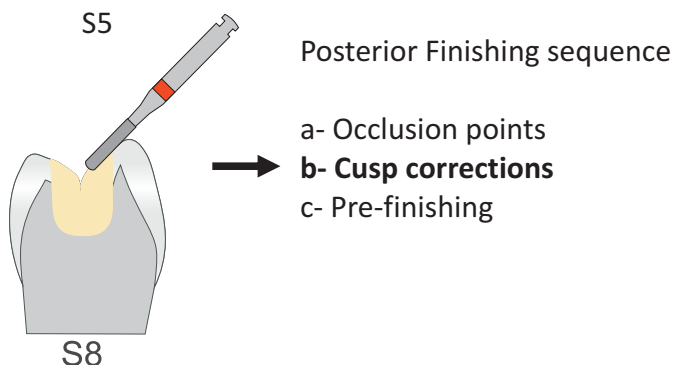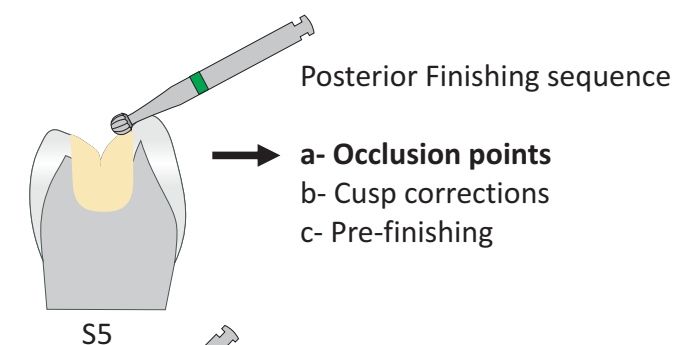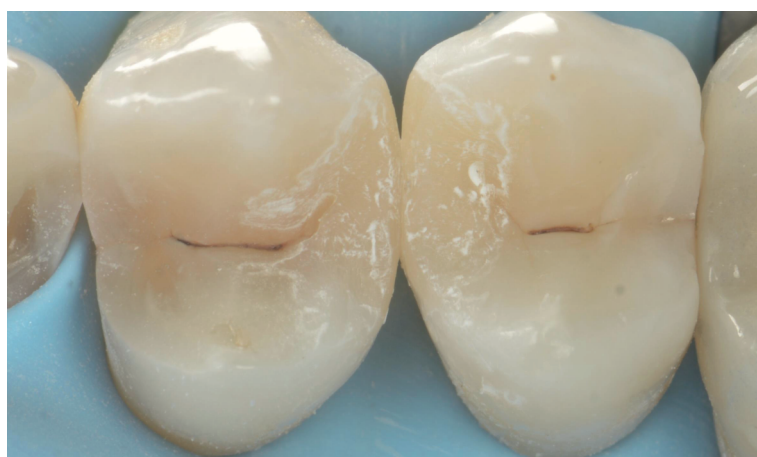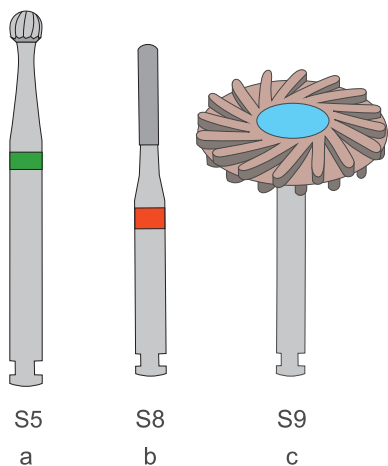

Styleitaliano

ABOUT STYLEITALIANO

Feasible, teachable and repeatable



The Style Italiano movement emerged from a simple idea of two friends – Walter Devoto and Angelo Putignano. The recognized experts in conservative and aesthetic dentistry had observed that materials allowing for the creation of beautiful direct restorations were available, but there was a lack of appropriate concepts for their successful use.


Anyone who wanted to imitate natural teeth in a perfect way needed artistic skills, creativity and lots of time to practice.

Since time is a scarce resource in every dental office, the two friends decided to develop simplified techniques that enable every dentist to obtain impressive results with direct composite materials.

Their practical suggestions are based on clear and precise ideas that help dentists face their own daily challenges in their dental office.



WALTER DEVOTO

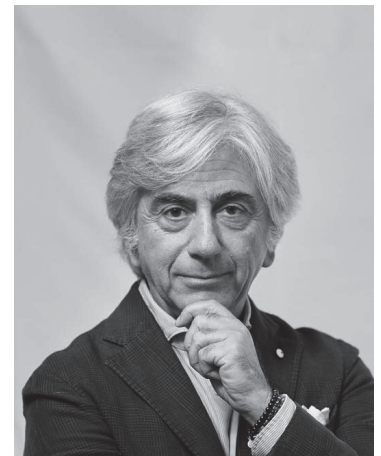

ANGELO PUTIGNANO



ABOUT HORICO

The company HOPf, Ringleb & CO. was founded by Ing. Paul Hopf in Berlin in 1918. Paul Hopf invented the first separating strips for the dental market. HORICO started its industrial production with discs and strip. Few years later rotary instruments were added to the range.

Famous inventions of HORICO:

- Steelcarbo®: Discs and strips
- Diaflex®: Flexible diamond disc
- Diatrepano®: Instrument for root canal opening
- Superdiaflex®: Diamond disc with a thickness of only 0.1 mm.

Today, HORICO is one of the world leaders in the manufacture and development of strips and rotary instruments (diamond and carbide) for dentistry and dental technology.

HORICO is a successful and innovative family business, run by the founder's great-granddaughters in the fourth generation and celebrating its 100th anniversary in 2018.



HOPF, RINGLEB & CO. GMBH & CIE.

Gardeschützenweg 82 - 12203 Berlin Germany
Tel.: +49(0)30 830 00 30 - Fax: +49(0)30 833 29 95
E-Mail: info@horico.de - Internet: www.horico.de



PR-KIT SI-DAY / 2026